

## BTS256-LED

<http://www.gigahertz-optik.de/en-us/product/BTS256-LED>

Product tags: VIS , Spectral Measurement



## Description

The photometric specifications of LEDs must commonly meet very high tolerance requirements even for non-specialist applications such as general and automotive lighting. This is often a problem since the manufacturing tolerances of LEDs can be higher than those permitted in the applications. The tolerance limits offered by LED manufacturers' intensity and color based binning are only applicable if the operating conditions are similar to those in the binning tests. Therefore, manufacturers incorporating LEDs into their products require devices that can accurately measure the precise in-situ photometric performance of LEDs.

### Compact spectral radiometer

The compact BTS256-LED enables you to conveniently measure the luminous flux, spectrum, color, and color rendering indices of single LEDs. One special feature is the conical measurement port of the device. The ability to perform measurements of onboard LEDs makes it possible to also include thermal effects in the measurement. The luminous flux, color, color rendering indices and spectrum of an LED are all typically measured within a few seconds. The device is therefore ideal for inspection of incoming products as well as the quality control in production processes. It can also be very useful in the design department.

The BTS256-LED comes in a compact aluminum housing and offers all functions that are necessary for precise measurement of the luminous flux, spectrum, color, and color rendering indices. For greatest accuracy and versatility, this device is based on a BiTec light sensor that consists of a V-lambda-filtered Si photodiode and a spectrometer unit that has a CMOS diode array. Si photodiodes are unsurpassed in terms of dynamic range, linearity, and speed. The CMOS diode array-based spectrometer guarantees precise measurement data of the luminous spectrum which is used to determine the color values. The combination of the two detectors enables mutual correction (see article on the [BTS technology](#)) for greater precision. This also makes it possible to perform accurate, time-synchronized measurements, e.g., of PWM signals. One cutting-edge feature of the BTS256-LED is its remote-controlled shutter for dark current compensation of the array as well as its software-controlled auxiliary lamp for compensation of light absorbed by the measurement samples (self-absorption correction). The remote control takes place via a USB 2.0 interface using the supplied [S-BTS256](#) software.

### Calibration


One essential quality of photometric devices is their precise and traceable calibration. Calibration of the BTS256-LED is performed in Gigahertz-Optik's calibration laboratory that is accredited by DAkkS (D-K-15047-01-00) for the *spectral responsivity* and *spectral irradiance* according to ISO/IEC 17025. The device has two calibrations: one is done using a specially developed reference lamp offering 2pi illumination which enables precise measurement of the luminous flux of diffusely emitting LEDs. The second calibration is for sources that have narrower illumination characteristics.

### Options for the BTS256-LED

- [Software development kit](#) to enable users to integrate the device in their own software
- Extension to the [BTS256-LED Plus Concept](#) (for illuminance and luminous intensity) using other components




*BTS256-LED for measurement of the luminous flux, spectrum, color, and color rendering index of single LEDs*

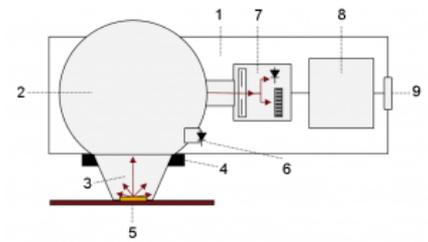


- 1) BTS256-LED housing
- 2) 50mm integrating sphere with synthetic coating
- 3) Conical measurement port
- 4) Precision bayonet Mount
- 5) Test LED on a circuit board (device under test)
- 6) Remote-controlled auxiliary lamp
- 7) BiTec sensor with Si photodiode, CMOS diode array spectrometer and shutter
- 8) Microprocessor
- 9) USB 2.0 interface




*The conical measurement port is placed over the test LED and detects all the radiation in a 2pi space*

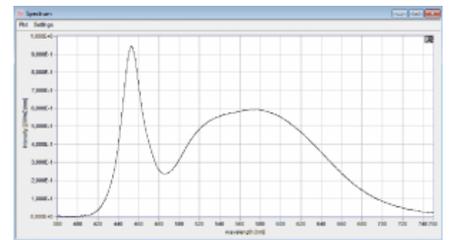


Graphical view of the spectrum

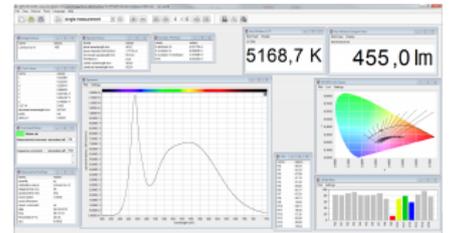

S-BTS256 user software with modular desktop setup

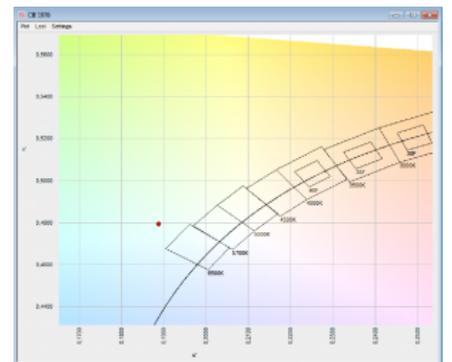

CIE 1976 chromaticity table with binning fields

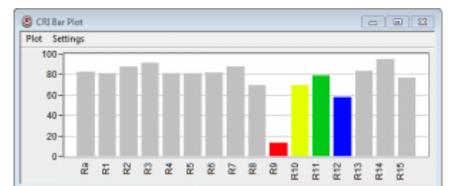

CRI Bar Plot

## Specifications

### General

Short description

Spectroradiometer for measurement of the luminous flux, spectrum, color, and color rendering indices of single LEDs

Main features	Compact measurement device with internal integrating sphere, BiTec light sensor, remote-controlled auxiliary lamp and shutter. Fast data logger for the luminous flux. Software
Measurement range	Luminous flux: 10 mlm - 1100 lm, spectral range: 360 nm - 830 nm, bandwidth: 5 nm with optical bandwidth correction according to CIE 214
typical applications	Goods-in inspection of individual LEDs, quality assurance of assembled LEDs in production processes, Research and development testing.
Calibration	For diffuse emitting and narrow beam LEDs. Factory calibration. Traceable to international calibration standards.
<b>Product</b>	
Calibration uncertainty	± 8 % for luminous flux
Sensor	Bi-Technology sensor with a photometric broadband detector and a array spectrometer. Integrated aperture for automatic dark signal adjustment.
Input optics	Integrating sphere with synthetic ODM98 coating and protective window at the sphere port. Cone adapter coated with ODP97 for radiation absorption. 10 mm diameter measurement port. LED auxiliary lamp.  Adapter change effect ± 0.5 % Max. xy responsivity deviation of the 10mm measurement port ± 2 % Max. z responsivity deviation of the 10mm measurement port ± 2 % (1 mm to 11mm)
<b>Spectral Detector</b>	
Chip	CMOS diode array
spectral range	(360 - 830) nm
Optical Bandwidth	5 nm
Data Resolution	1 nm
Integration Time	(5.2 - 30000) ms
Shutter	Automatic aperture for dark signal measurements with the same integration time as that of light measurements. Aperture delay = 100ms .
typical measurement time	1100 lm ≤ 5ms (white light)  10 mlm ≤ 30s (white light)
Peak wavelength	± 0.5 nm
Dominant wavelength	± 1 nm
	± 0.0001 (Standard illuminant type A) ± 0.0002 (LED)
	± 0.002 (Standard illuminant type A) ± 0.005 (typ. LED)
CCT Measurement range	(1700 - 17000) K
ΔCCT	± 50K (standard illuminant type A) ± 3% (depending on the LED spectrum)
CRI (color rendering index)	Ra and R1 to R15
Stray Light	6E-4 (Blue LED)  6E-4 (Green LED)  6E-4 (Red LED)


1E-3 (White LED)

### Integral Detector

max. luminous flux	70000 lm
Filter	Spectral responsivity with fine CIE photometric matching. Online correction of the photometric matching through spectral measurement data (spectral mismatch factor correction).
f1'	≤ 6 % (uncorrected) ≤ 3 % (f1' a*(s <sub>z</sub> (λ)) respectively F*(s <sub>z</sub> (λ)) corrected by spectral data, done automatically by BTS technology)
ADC	12Bit
Measurement time	(0.1 - 6000) ms
Noise equivalent luminous flux	0.05 mIm

### Graphs

spectral responsivity

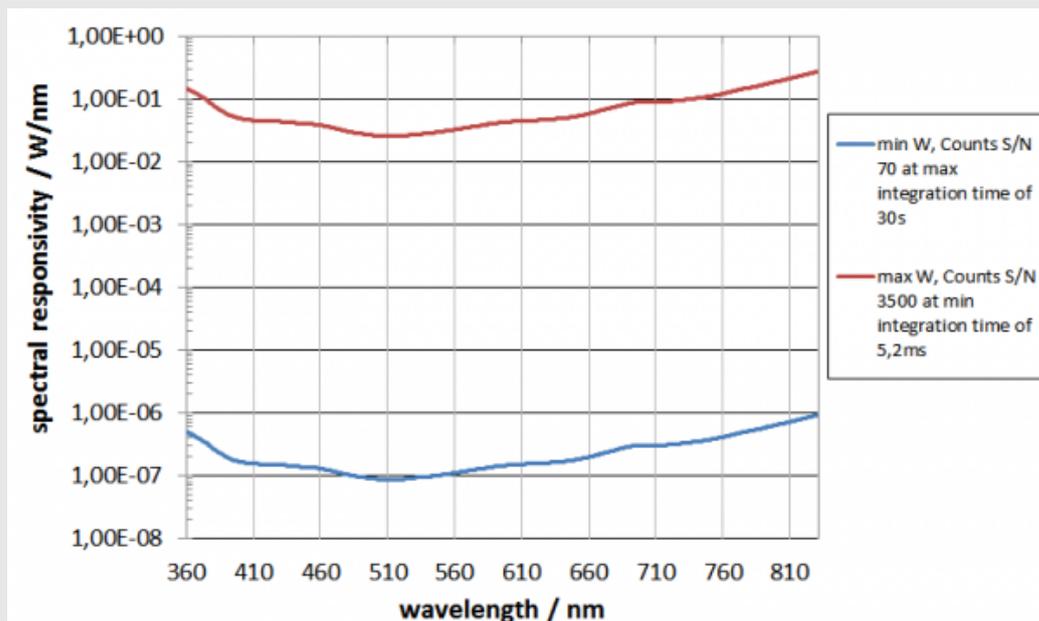

### Miscellaneous

Microprocessor	16Bit, 25ns instruction cycle time
Power Supply	5VDC to 7VDC, 250mA peak during capacitor charging of the auxiliary lamp
Interface	USB 2.0 (Type B USB port)
temperature range	Operation: (10 to 30) °C Storage: (-10 to 50) °C
Dimensions	160 mm x 85 mm x 60 mm (Length x Width x Height)
Weight	500 g
Transport case	Plastic hard-top casing, 333 mm x 280 mm x 70 mm, 650g

#### Option: 210mm Integrating Sphere (UMBB-210)

Luminous flux measurement range (integral measurement)	(0.35 - 35000) lm
Sphere diameter	210 mm
Calibration	Luminous flux: ± 8%

#### Option: 500mm Integrating Sphere (UMBB-500)

Luminous flux measurement range (integral measurement)	(1.2 - 120000) lm
Sphere diameter	500 mm
Calibration	Luminous flux: $\pm 8\%$
<b>Option: 1000mm Integrating Sphere (UMTB-1000-HFT)</b>	
Luminous flux measurement range (integral measurement)	(4 - 400000) lm
Sphere diameter	1000 mm
Calibration	Luminous flux: $\pm 8\%$
<b>Option: Irradiance (DA)</b>	
Measurement range	Illuminance: (0.2 - 25000) lx
Calibration	$\pm 4\%$
<b>Option: Goniometer (GB-GD-360-RB40)</b>	
Luminous intensity measurement range (integral measurement)	(2E-1 - 2E8) cd ; by 1m measurement distance
Calibration	Luminous intensity: $\pm 4\%$

## Configurable with

Produktname	Product Image	Description	Show product
BTS256-LED-IB		<p>Compact Bi-Tec measurement device for the measurement of ILED-B.</p> <p>Features: Bajonett adapter with ILED-B geometry according CIE 127, spectral data, color temperature, CRI, chromaticity coordinates, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-IB">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-IB</a>
BTS256-LED-DA		<p>Compact Bi-Tec measurement device for the measurement of illuminance and luminous flux.</p> <p>Features: Bajonett adapter with diffusor for the BTS256-LED, +/- 30° cosine corrected field of view, spectral radiant power, color temperature, CRI, chromaticity coordinates, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-DA">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-DA</a>
ISD-21-BTS256-LED		<p>System for the luminous flux and light color measurement of individual 2π light emitting diodes up to 63.5mm.</p> <p>Features: Integrating sphere with 21cm diameter and auxiliary lamp, compact spectral light meter with Bi-Tec sensor for accurate measurement of the luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-21">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-21</a>
ISD-50-BTS256-LED		<p>System for the luminous flux and light color measurement of individual 2π light emitting diodes up to 100mm.</p> <p>Features: Integrating sphere with 50cm diameter and auxiliary lamp, compact spectral light meter with Bi-Tec sensor for accurate measurement of the luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-50-V01">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-50-V01</a>

Produktname	Product Image	Description	Show product
ISD-50HF-BTS256-LED		<p>System for the luminous flux measurement of <math>4\pi</math> (optional <math>2\pi</math> and <math>4\pi</math>) light fixtures inside a sphere.</p> <p>Features: Integrating sphere with 50cm hemispherical shell for opening and closing and auxiliary lamps. Light meter with Bi-Technology sensor for accurate measurement of luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/ISD-50HF-BTS256-LED">http://www.gigahertz-optik.de/en-us/product/ISD-50HF-BTS256-LED</a>
ISD-100HF-BTS256-LED		<p>System for the luminous flux measurement of <math>2\pi</math> and <math>4\pi</math> light fixtures inside a sphere.</p> <p>Features: Integrating sphere with a 1000 mm diameter, hemispherical shell for opening and closing, auxiliary lamp. Light meter for accurate measurement of luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, User software, etc. Options: quadripole test sockets.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-100HF-V01">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-100HF-V01</a>
ISD-100HF-V02-BTS256-LED		<p>System for the luminous flux measurement of <math>2\pi</math> and <math>4\pi</math> light fixtures inside a sphere.</p> <p>Features: Integrating sphere with a 1000 mm diameter, extra measurement ports for <math>2\pi</math> luminaires with diameters of up to 254mm and auxiliary lamp. Light meter for accurate measurement of luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, User software, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-100HF-V02">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-ISD-100HF-V02</a>
GB-GD-360-RB40-2-BT S256-LED		<p>Goniometersystem for the luminous flux and light color measurement of individual <math>2\pi</math> light emitting diodes with an adjustable measurement distance up to 2000mm.</p> <p>Features: Goniometer, compact spectral light meter with Bi-Tec sensor for accurate measurement of the luminous flux, spectral radiant power, CCT, CRI, chromaticity coordinates, IES, EULUMDAT, etc.</p>	<a href="http://www.gigahertz-optik.de/en-us/product/BTS256-LED-GB-GD-360-RB40-2">http://www.gigahertz-optik.de/en-us/product/BTS256-LED-GB-GD-360-RB40-2</a>
S-SDK-BTS256		Software Development Kit for BTS256 variants.	<a href="http://www.gigahertz-optik.de/en-us/product/S-SDK-BTS256">http://www.gigahertz-optik.de/en-us/product/S-SDK-BTS256</a>

## Purchasing information

Article-Nr	Modell	Description
<b>Product</b>		
15308420	BTS256-LED	Measurement device, BTS256-LED-CA10 cone adapter, USB cable, hard-top casing, operation manual, S-BTS256 software, calibration certificate.
<b>Re-calibration</b>		
15300226	K-BTS256-LED-I	Recalibration of the BTS256-LED Tester. Only possible with the 10mm cone adapter
<b>Software</b>		
15298218	S-SDK-BTS256	Software Development Kit for the implementation of the BTS256 or variants into custom made software