

LS-OK30-FH


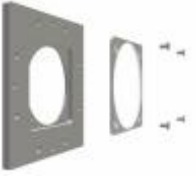

<http://www.gigahertz-optik.de/en-us/product/LS-OK30-FH>

Product tags:

Description

The LS-OK30-FH enables installation of neutral and/or correction filters between the LS-OK30 light source and the integrating sphere. This makes intensity change in different levels or spectral influence on the luminous spectrum possible. Due to the extremely high temperatures of the halogen lamps, not all filter designs are appropriate. The 50mm diameter filter with a 3mm maximum width is placed in type LS-OK30-FHI filter pocket. This enables tool free filter change. The filter holder can be combined with the variable baffles.

Configurable with

Produktname	Product Image	Description	Show product
LS-OK30		Light source for use with the UM series modular construction integrating spheres. Features: External source for LH-F series halogen lamps up to 100W. Diffuse reflector with synthetic coating. Fan for active air-flow lamp socket cooling.	http://www.gigahertz-optik.de/en-us/product/LS-OK30
UMPF-LSOK30		Port frame for use with the lamp LS-OK30. Features: Optional diffuse window for use between LS-OK30 and the sphere.	http://www.gigahertz-optik.de/en-us/product/UMPF-LSOK30
LS-OK30-VA		Variable aperture for use with lamp LS-OK30. Features: High resolution manual drive.	http://www.gigahertz-optik.de/en-us/product/LS-OK30-VA

Purchasing information

Article-Nr	Modell	Description
Product		
101011	LS-OK30-FH	Filter holder for LS-OK30-FHI filter frame inserts.
Accessories		
101014	LS-OK30-FHI	Filter frame insert for LS-OK30-FH.